



MiniMe[®] 2

Paediatric nasal NIV mask



Critical Care ▪ Patient Interfaces



MiniMe®2

Featuring AIR^ogel® technology

The MiniMe®2 mask is a single-patient use nasal NIV mask for use in the hospital or homecare environment. MiniMe®2 is available in two sizes, both with vented and non-vented options.

Designed for children ages 2-12, the MiniMe®2 mask features AIR^ogel® with Advanced Cushion Technology™, providing your smallest patients the ultimate fit and maximum comfort.

Custom fit technology

The shell and cushion of the mask can easily be moulded to fit individual needs. Once fixed it will hold its position



AIR^ogel®

The AIR^ogel® makes the MiniMe®2 a gentle fitting mask

Touchless Spacebar™

Unique design eliminates any contact with the patient's forehead

Size indicator

The masks are available in two sizes, small/medium and medium/large

Advanced Cushion Technology™


The ergonomic cushion technology delivers an effective seal for the most comfortable fit

Vented and Non-vented versions available

For CPAP and bilevel therapy in children

360° rotating elbow

For optimal tube positioning and easier mobility

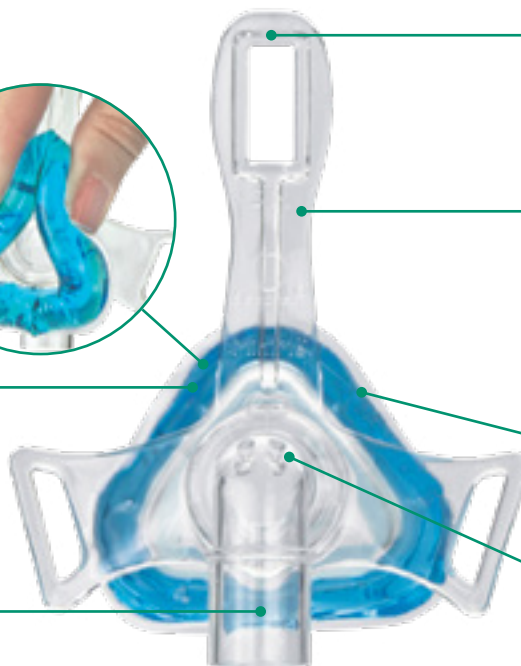

Code	Description	Colour coded elbow	Box. Qty
The masks below should be used with single limb breathing systems that do not contain a CO ₂ leak port			
2370000	MiniMe®2 Paediatric nasal vented NIV mask small/medium	○	10
2371000	MiniMe®2 Paediatric nasal vented NIV mask medium/large	○	10
The masks below should be used with single limb breathing systems that have a CO ₂ leak port			
2270000	MiniMe®2 Paediatric nasal non-vented NIV mask small/medium	●	10
2271000	MiniMe®2 Paediatric nasal non-vented NIV mask medium/large	●	10

Replacement headgear for small/medium MiniMe®2

Code	Description	Box. Qty
2274000	55048 Headgear	20

Replacement headgear for medium/large MiniMe®2

Code	Description	Box. Qty
2275000	EZ Fit Headgear, packaged	20

Make an enquiry

MiniMe®2 and AIR^ogel® are registered trademarks of Sleepnet Corporation.

Touchless Spacebar and Advanced Cushion Technology are trademarks of Sleepnet Corporation.

Please think before you print

Save energy and paper. If you must print this information sheet please print it double sided.

IS1.49 ▪ Issue 1 02.19



Intersurgical Ltd, Crane House, Molly Millars Lane, Wokingham, Berkshire, RG41 2RZ, UK
T: +44 (0)118 965 6300 F: +44 (0)118 965 6356 info@intersurgical.com www.intersurgical.com

UK ▪ Ireland ▪ France ▪ Germany ▪ Spain ▪ Portugal ▪ Italy ▪ Benelux ▪ Sweden ▪ Denmark ▪ Lithuania ▪ Russia ▪ Czech Republic
Turkey ▪ South Africa ▪ China ▪ Japan ▪ Taiwan ▪ Philippines ▪ USA ▪ Canada ▪ Colombia ▪ Australia