



Neonatal Resuscitation

with variable PEEP



Critical Care ▪ Breathing Systems for Critical Care

Neonatal resuscitation breathing systems with variable PEEP

A range of 10mm Neonatal resuscitation products which incorporate an adjustable orifice control on the patient elbow to provide variable PEEP.

Flow through the breathing systems is set in line with the patient's inspiratory peak flow requirement. Peak inspiratory pressure is also set on the resuscitation device by the clinician, so PEEP can now be adjusted adjacent to the patient using the click adjustment on the PEEP cap if required.

The versatile 22M/15F patient connector not only fits all interfaces but swivels to aid patient comfort and increase seal efficiency.

These breathing systems can be used in NICU and obstetric units with the following resuscitation platforms and gas drivers:

- Dräger®
- Hilrom Airshields
- GE® Datex-Ohmeda®, Giraffe®, Panda®
- Viamed Tom Thumb®
- Fisher & Paykel® Neopuff®


The breathing systems are available in a variety of lengths to accommodate various types of equipment and positioning of the patient.

0.8m, 1.2m and 2m tubing options

Offering more choice for different ventilatory systems

Click adjustment on PEEP cap

No accidental adjustment during operation

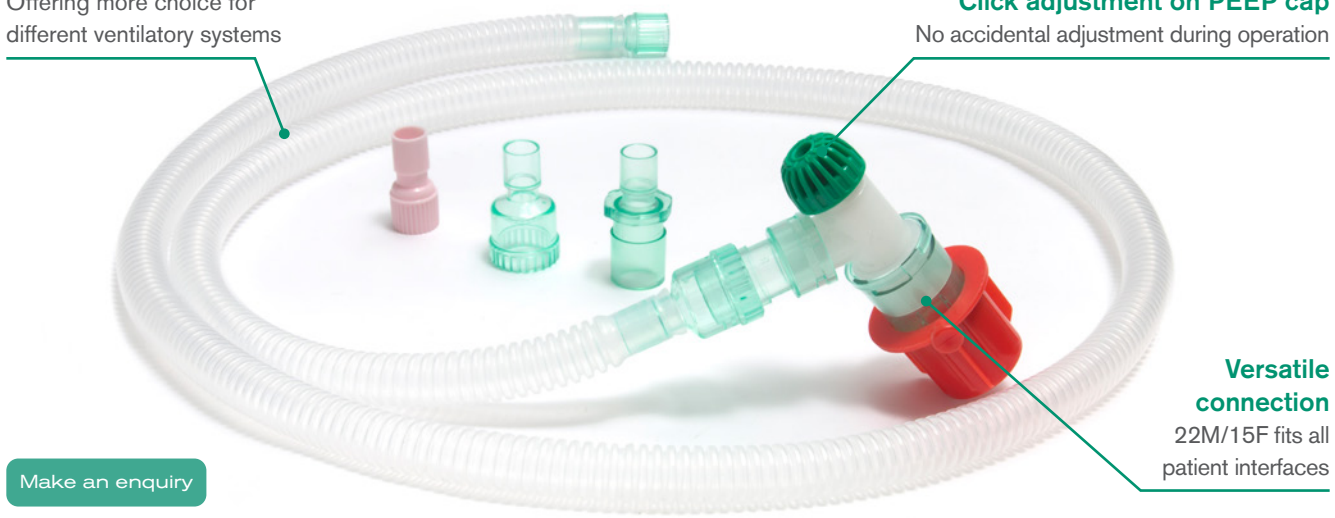

Make an enquiry



Code	Description	Box Qty.
6430000	10mm Flextube neonatal resuscitation breathing system with variable PEEP, double swivel elbow, 0.8m	15
6431000	10mm Flextube neonatal resuscitation breathing system with variable PEEP, double swivel elbow, 1.2m	15
6432000	10mm Flextube neonatal resuscitation breathing system with variable PEEP, double swivel elbow for Neopuff®, 1.2m	15
6432001	10mm Flextube neonatal resuscitation breathing system with variable PEEP, double swivel elbow for Neopuff®, 2m	10
6433000	10mm Flextube neonatal resuscitation breathing system with variable PEEP, double swivel elbow, Neopuff® and universal connectors, 1.2m	15

Dräger is a registered trademark of Drägerwerk AG & Co. KGaA. GE is a registered trademark of General Electric Company. Datex-Ohmeda, Inc., is a General Electric company. Datex-Ohmeda is a registered trademark of GE Healthcare Finland Oy. Giraffe and Panda are registered trademarks of Datex-Ohmeda, Inc. Tom Thumb is a registered trademark of Viamed, Ltd. Fisher & Paykel is a registered trademark of Fisher & Paykel Healthcare Limited. Neopuff is a registered trademark of Fisher & Paykel Healthcare Limited.



Intersurgical Ltd, Crane House, Molly Millars Lane, Wokingham, Berkshire, RG41 2RZ, UK
 T: +44 (0)118 965 6300 F: +44 (0)118 965 6356 info@intersurgical.com www.intersurgical.com



The manufacturer Intersurgical Ltd is certified to ISO 9001:2015, ISO 13485:2016 and ISO 14001:2015

Please think before you print
 Save energy and paper.
 If you must print this information sheet please print it double sided.